

Before Use



When using for the first time, press the silicone valve vertically to open the slits before suctioning.



When suctioning, keep the hand pump and the valve vertical.



The extension tube will increase the range of self-application when suctioning.



It is recommended to warm up the receiver's body before cupping.



Applying skin moisturizing lotion or oil will help suction and adhesion.



Made in Korea for Meridius Medical GmbH

EU Representative : Meridius Medical Europe Ltd.
Unit 3D, North Point House, North Point Business Park, New Mallow Road, Cork, T23 AT2P, Ireland

Importer: Meridius Medical GmbH
Heiderhöfen 23, 46049 Oberhausen, Germany
Tel: +49 (0)208 20580580 eu@meridiusmedical.com



User's manual Meridius Cupping Set



Move to
Homepage

ME30

contains 20 square cups
and 10 curved cups.



Components



1 Hand Pump



20 square cups



2 CO1 cups



2 CO2 cups



1 Extension Tube



2 CO3 cups



2 CO4 cups



RO4 cup TR cup


Precautions for use

Please read User's manual carefully before use.

01. For pregnant person, do not apply directly to the abdomen. For patients with low blood pressure or weakness, the duration and intensity of use must be increased step by step.
02. Apply skin moisturizing lotion or oil to protect the skin and sustain vacuum.
03. If patient feel dizzy, experience severe pain, or show a blister, stop immediately and treat according to safety code.
04. Do not apply on eyes and ears.

Not a faulty valve

The silicone valve has three slits, through which air travels. Before first time use or if unused for a few weeks, press the silicone valve vertically to open the slits.

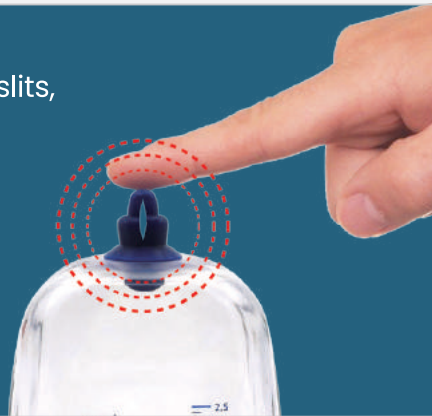

Hypothetical suggestion to optimize stimuli with NEW rounded square cups.

01. Attach the primary cup to the place that requires optimum stimuli.
02. Attach the auxiliary cups in 4 or 8 directions around the primary cup marked in red in the pictures.

Relatively evenly applied multiple square cups are asserted to optimize negative pressure throughout wide surface, therefore magnitude of stimuli and clinical effects.



Versatile cups to fit various body contours (n=10)

Five different shapes and sizes of versatile cups are to be applied to the neck, elbows, wrist, knees and ankles.



Cups dry well even when stacked

The square cups are designed to allow drainage and ventilation even when stacked.



Detachable and replaceable silicone valve

Detachable and replaceable silicone valve allows convenient and efficient disinfection.



Disinfect the cups according to local standard

Please sterilize them using cotton wool balls soaked with methylated spirits. It also can be used in Autoclave. Make sure to remove the silicone valves.



Enhanced suction with a double ring design

Our double ring design inside the pump piston enhances the suction power.

